

# Salsa & Salsa M

## Options and Accessories





# LIVE WITHOUT LIMITS



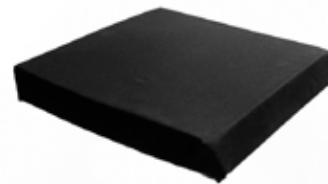
SEATING.....	4
JAY COMFORT SEATING.....	5
JAY CLINICAL CUSHIONS AND BACKS.....	9
BACKREST RECLINE OPTIONS.....	11
POWER SEATING MODULES.....	12
ARMRESTS .....	13
LEGRESTS.....	13
Swing Away & Elevating Legrests.....	13
Footplates.....	15
Legrest accessories.....	15
HEADREST.....	16
CONTROL SYSTEMS.....	17
VR2 Controls.....	17
Control Boxes.....	17
R-NET Controls.....	18
R-NET optional controls/electronics.....	20
Lights & Indicators.....	21
WHEELS .....	22
BATTERIES & CHARGERS .....	23
FRAME COLOUR .....	24
ACCESSORIES .....	25
QUICKIE WEB IQ ACCESSORIES .....	27

# Seat



Firm seat board

SAL040207



5 cm (2") standard foam cushion

41 cm wide SAL040200

46 cm wide SAL 040201



Tension adj. sling backrest

SAL050001

Backrest height:

SAL050002 - 46 cm

SAL050003 - 48.5 cm

SAL050004 - 51 cm



Standard comfort seat cushion &  
Backrest

41 cm width

46 cm width

51 cm width

SAL040208



# JAY Comfort Seating

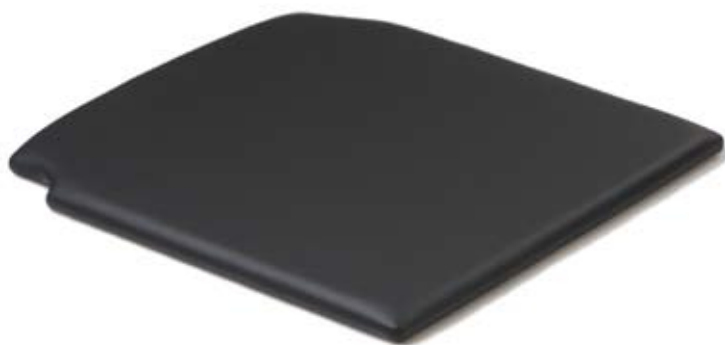
The Jay Comfort Seating system provides clients with stylish seating that has been designed to maximise comfort whilst meeting mild to moderate clinical requirements.

## Jay comfort seat pan

A firm padded board is a solid base for any seat cushion.

This seat pan with easy clean upholstery is available in different sizes

SAL040217



---

## Jay comfort seat pan with 5 cm (2") threshold

The comfort seat pan with added thresholds is designed for clients with their own cushion in mind.

The added thresholds give the benefit of stopping the movement between the cushion and the seat pan, giving a stable base for clients.

SAL040218





# JAY Comfort Seating

## Jay comfort standard cushion - stretched rubber upholstery

This cushion has a 12 mm ( 1/2") beech wood base with a slightly contoured foam cushion, featuring a 2 cm (3/4") seat angle.

SAL040219



---

## Jay optimum comfort - stretched rubber upholstery

The Optimum comfort cushion includes a 12 mm (1/2") beech wood base with contoured foam cushion.

SAL040220





### Jay standard velour backrest - hard back

This anatomically shaped backrest features a luxury velour finish for added warmth and comfort and an easy clean rubber upholstery to the sides and rear.

SAL050010



Rear view

### Jay standard velour backrest - vertical elastic straps

This Standard Velour comfort backrest features flexible elastic straps instead of the hard wooden structure featured in the hard back version shown above.

The flexible elastic straps provide a softer and more comfortable ride whilst still providing adequate support.

SAL050011



Rear view

# JAY Comfort Seating

## Jay standard velour backrest - horizontal adjustable straps

This standard velour comfort backrest features tension adjustable straps which can be adjusted and shaped to the client's back, providing excellent support and flexibility.

SAL050012



Rear view



## Jay deep contour backrest - active

This comfort backrest features a deeper contour, anatomically shaped to provide support and stability around the lower back and is cut out around the shoulders so as not to restrict active clients.

Tension adjustable straps also provide additional support and flexibility for the client.

SAL050014



Rear view



## Jay deep contour backrest - passive

This comfort backrest features a deeper contour, anatomically shaped to provide support and stability across the entire back for clients with restricted movement.

Tension adjustable straps also provide additional support and flexibility for the client.

SAL050013



Rear view







# Jay clinical cushions and backs

JAY® offers a complete range of cushions and backs to meet individual user needs. The comfort cushions, for users with moderate to low skin protection and positioning needs. Skin protection cushions for moderate skin protection and positioning needs and the Positioning & Skin protection cushions, for user with more complex positioning needs and or high skin breakdown risk.

For more information please refer to the JAY® cushions Brochure and Order form



Basic



Soft Combi P



Easy Visco



Xtreme



Easy Fluid



Active



J3



J2



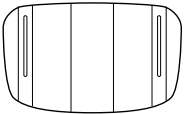




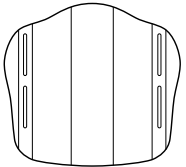




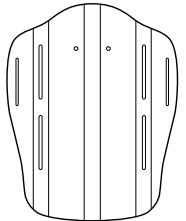




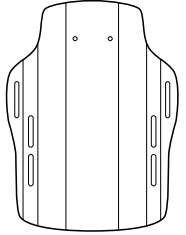






GS

# Jay Backs

JAY® backs are designed to provide posterior lateral pelvic stabilization along with enhanced trunk stability. They are lightweight and adjustable, allowing for easy installation. JAY® backs reduce fatigue and maximize function to provide clients the comfort and stability to go anywhere with confidence.

For more information please refer to the JAY® backs Brochure and Order form

 The contour depths and back heights are available in the following combinations					
 Lower thoracic (LT)		Shallow Contour (SC)	Mid Contour (MC)	Deep Contour (DC)	Posterior Deep Contour (PDC)
					
Headrests available for SC and MC					
 Mid thoracic (MT)					
Lateral, headrest and harnessis available for SC, MC, DC and PDC					
 Upper thoracic (UT)					
Lateral, headrest and harnessis available for SC, MC, DC and PDC					
 Shoulder height (SH)					
Width available:		31, 36, 41, 46 and 51 cm	31, 36, 41, 46 and 51 cm	36, 41, 46 and 51 cm	36, 41, 46 and 51 cm

# Backrest Recline Options

The Quickie Salsa & Salsa M has a recline option to suit everyone. Via tool adjustment, carer adjustment or client adjustment and via a manual or power solution.



## Mechanical fixed recline backrest

The standard rehab manual adjustable back provides a fixed back angle from  $-4^{\circ}$  to  $16^{\circ}$ .

Adjustment is still quick and easy. Simply remove the fixing bolts, set the desired angle and refit the bolts.

SAL050100



## Power recline

Ideally suited for clients who want to adjust the back angle with no assistance from a carer.

Powered from the hand control, the back can be reclined from  $-5^{\circ}$  to  $30^{\circ}$ .

SAL050102

---

## Manual recline - gas spring system

If your client requires a back that can be reclined whilst in the chair, then the gas spring system is for them.

With the client seated in their SALSA, carers can recline this back from  $-5^{\circ}$  to  $30^{\circ}$  with the press of a lever.

SAL050101



# Seat to base interface

## Power Seating Module

The Salsa & Salsa M power chair offers simple and versatile solutions for positioning. A range of options are available. Manual preset tilt of 0 - 9 degrees, power tilt, power lift or a combination of both powered lift and tilt.

Simple modular units can be added or taken away as desired. This capability ensures that the Salsa & Salsa M power chair is future-proof.

### Power tilt module (0°- 30°)

The Quickie Salsa power tilt provides clients with 0° to 30° of tilt adjustment controlled through most input devices.

SAL040212

### Power lift & tilt combi module

The Quickie Salsa lift & tilt module provides clients with 12" (30cm) of lift and 0° to 30° of tilt adjustment, powered through most input devices.

SAL040214



Manual tilt (3° preset tilt)

SAL040211



Power tilt module (-5° - 25°)

The Quickie Salsa powered tilt -5° to 25°

SAL040213



Protective cover for powered tilt module

SAL040215

Protective cover for powered lift / tilt combi module

SAL 040216



## Armrests

The Quickie Salsa & Salsa M offers a comprehensive range of armrests and accompanying sideguards to suit everyone. Once you have chosen your armrest we have a wide range of armpads from standard desk length pads to arm and hand pads designed specifically for Hemiplegic clients.



Reclining armrest function

SAL 060001

Flip up, height and  
angle adjustable  
c/w 385 x 60 mm pada

SAL060000



## Legrests

Articulating, elevating and different angles of hangers, together with a wide range of footplates to suit every client preference.



70° hanger

### Swing-away hangers

Push the button to swing away and detach. The range of swing away legrests is available in all frame colours.

SAL060200



### Extension tube for s/a legrests:

lower leg length : 34- 48 cm  
( 13"-19")

SAL060203



### Auto - Connectors for Powered s/a legrests

Fixed width legrest receiver  
The standard legrest receiver  
offering no width adjustment.

SAL060208



### Manual elevating legrest

The manual elevating legrest can be easily adjusted by the flip of a lever.

Swing-away, detachable and complete with knee pads as standard.

SAL060201



### Powered elevating legrest

The powered elevating legrest is powered from the hand control and gives even the most challenged client the ability to adjust the legrest.

Swing-away, detachable and complete with knee pads as standard.

SAL060208 Auto - Connectors available

SAL060202



### Powered articulating legrest

Ideal for clients with severe clinical requirements and powered from the hand control, the articulating legrest automatically extends as it elevates.

SAL060208 Auto - Connectors available

SAL060205



### Fixed centre mount legrest

Adult: Lower Leg Length 37 cm to 42 cm  
Kids : Lower Leg Length 20.4 cm - 31.6 cm

Only available for Salsa M

SAL060204

SAL 060206



## Footplates



### Aluminium angle adjustable

Made from aluminium this footplate will benefit clients looking for a more stable platform, and the ability to adjust the angle of the footplate.

SAL060300



### Footboard

The one piece, angle adjustable footboard is constructed from lightweight aluminium giving clients the most stable platform available, whilst still retaining the benefit of detaching one side and flipping it up.

SAL060301

## Legrest accessories



### Legrest strap



### Heel loop

Standard on the composite flip up, angle adjustable and aluminium flip up footplates



### Pro seating comfort calf pad

The comfort seating calf pad is made from a thick extra soft foam and an easy to clean cover for added comfort and convenience.

SAL060207

### Extension brackets for std. calf pads

SAL060209

# Headrest

## Comfort headrest

With adjustable side pads ,can be adjusted to give optimum support to clients

SAL090002



## Comfort Seating - standard headrest

This Comfort seating headrest provides a soft velour finish for warmth and easy-clean rubber upholstery.

SAL090014



## Comfort Seating - flexible headrest

This comfort headrest provides the same benefits as the standard headrest with the added ability to adjust the flexible wings around the client's head for added support

SAL090015



## Curved headrest 23 cm x 11 cm

SAL090025





## VR2 controls

If you just need a high performance base with a basic controller function the VR2 standard controller is your choice. It provides 90A motor output for enough torque to manage smooth or even rough terrains.

VR2 JSM with  
actuator control

SAL030050



VR2 JSM with  
actuator and light  
control

SAL030051



VR2 JSM  
Attendant only

SAL090016



VR2 4 actuator  
control box

SAL030054



VR2 dual control

SAL090004

## Control boxes

VR2 Lighting module

SAL030061



VR2 70A motor controller

SAL030063

VR2 90A Heavy duty motor controller

SAL030064



# PG R-NET controls



The PG R-Net control system breaks the mould with design, modularity and advanced performance, offering more functionality than any other control system on the market.

For the first time, dealers and therapists have the ability to programme and tailor drive and seating profiles to client needs, and assign multiple functions directly to hand control buttons making life easier for the user.

## R-Net LED Joystick, with Actuator Controls

Multiple profile hand control used for drive and actuator control.  
Assignable Button function built into this joystick.

SAL030057



---

## R-Net LED Joystick, with Actuator Controls & Lights

Multiple profile hand control used for drive, actuators and lighting control.

Assignable Button function built into this joystick.

SAL030058





## R-Net LCD Joystick

This joystick has an LCD screen providing clients with visual feedback. The colour display shows large icons in combination with programmable text (profiles & Modes)

With up to 8 programmable profiles and speed button to finely adjust speeds within profiles this joystick also has a clock, Speedometer, and odometer.

Complete with 2 separate 3,5mm Jack sockets for external switches (On/Off and Mode - Using a splitter and the PC programmer you can assign to other functions (4 in total)

Assignable Button function built into this joystick.



SAL090017

---

## R-Net Attendant only

SAL090018

## R-Net Dual Control

Compact multiple profile hand control mounted at the rear for an attendant.

Control profiles and seating functions with this joystick.

SAL090004



# R-NET optional controls / electronics

## R-Net Lighting Module (ISM-L)

This module is required when Lights & Indicators are desired.

SAL030059

## R-Net Seating & Lighting Module (ISM-6L)

This seating and lighting module is required to upgrade the system from one to a maximum of six powered seating functions, whilst also providing lights and indicators function.

SAL030060



---

## Bluetooth USB mouse receiver

SAL 090012



---

## R-Net Bluetooth PC Control Module

The bluetooth module provides users with the ability to control there PC mouse from the joystick or other input device. e.g. Chin Control

The users PC or laptop requires bluetooth functionality or a bluetooth receiver is required (RNT090506)

SAL090011







## Lights

Choose from our standard STVZO approved lights and indicators or exclusive to the Quickie Salsa & Salsa M, LED lights and indicators. These project a more focused and intense beam providing greater visibility. The LED Lights also draw less power from the battery and have a considerably longer life than standard lights.

### Standard lights

Front & rear lights and indicators

SAL020000



---

### Led lights

LED front & rear lights and indicators

SAL020001



# Drive wheels

The Quickie Salsa & Salsa M offers an extensive range of drive wheels and castors to suit most types of drive terrain and environment. Choose from a variety of tyre treads and sizes, available in solid, pneumatic or OKO filled (puncture resistant).

## Salsa

- SAL070000 8.5" solid castor , 12.5" Pneumatic Drive wheels
- SAL070001 8.5" solid castor , 12.5" Solid Drive wheels
- SAL070002 9" solid castor , 14" Pneumatic Drive wheels
- SAL070003 9" solid castor , 14" Solid Drive wheels
- SAL070004 9" solid castor , 14" Solid Knobbly Drive wheels
- SAL070005 9" solid castor , 14" Pneumatic Knobbly Drive wheels



Solid / Pneumatic Knobbly Drive wheels



Solid / Pneumatic Drive wheels

- SAL070010 Emergency Free wheel mechanism, with 8.5" solid castor, 12.5" solid Drive wheels
- SAL070011 Emergency Free wheel mechanism, with 9" solid castor, 13" solid Drive wheels

---

## Salsa M

- 7" solid front castor , 6" solid rear castor, 13" Pneumatic Drive wheels
- 7" solid front castor , 6" solid rear castor, 13" Solid Drive wheels



- SAL070011 Emergency Free wheel mechanism, 13" solid Drive wheels



## Batteries & Charges

The Quickie Salsa & Salsa M offers a choice of batteries and chargers.

Clients who are likely to use the Salsa & Salsa M for drive only and basic functions will find the 60Ah batteries suitable, but if the Groove is equipped with features like an enhanced user display, lift or tilt and lights and indicators

### Batteries

60 Ah maintenance free Only available in Salsa R

SAL030100

50 Ah gel batteries Only available in Salsa M

SAL030102

---

### Charges

220V / 8A automatic charger

SAL090200



# Frame colour



Salsa Colour Silver

SAL090300

---

Salsa M Colour Red

SAL090301







## Accessories

The Quickie Salsa & Salsa M offers an extensive range of accessories that have been designed to assist users with their day to day living.

Transporting tie down points



Rear view mirror

SAL090006



Crutch holder

SAL090003



5 cm (2") positioning strap

SAL0090005



# Accessories

## Swing-away tray table

The swing-away tray table provides clients with an easy to clean and stable surface to work or eat from. Once finished simply swing out of the way or detach.

SAL 090001 ( 41 - 48 cm)

SAL 090013 (46 -51 cm)



---

## Swing-away tray-mounted centre control

The swing-away tray-mounted centre control enables the hand control to be flipped over providing a flat surface when the tray is required and back around when the user requires access to the hand control.

SAL030062



---

## Kerb climber

Only available in Salsa R

SAL030401





## Quickie Web IQ Accessories

For the first time Quickie dealers are able to easily reconfigure and personalise the QuickieSalsa & Salsa M using the revolutionary Quickie IQ software kit and hand held programmer.

DTT with R-Net Programming cable

SAL090021

DTT with VR2 - Programming cable

SAL090020



DTT

### R-Net PC Programmer

The R-Net PC Programmer enables the dealer to re configure and programme a client's Quickie Groove allowing drive parameters and settings to be adjusted to provide the perfect drive tailored to that client.

Also provides diagnostics and test functions for service engineers.



SAL090010

### P&G PC Mobility programmer ( VR2 programming)

SAL090007



Web IQ is the world's first interactive online programming tool, enabling dealers to create, review and store unique customer profiles in a safe and secure environment.

Features include:

- Ability to download chair configuration and reconfigure the Power wheelchair - Review chair history
- Save drive profiles on a secure server (located at Sunrise Medical)
- Download files to a chair using Quickie Web IQ.

[www.sunmedwebiq.com](http://www.sunmedwebiq.com)



#### **Sunrise Medical GmbH & Co. KG**

Sunrise Medical S.L  
Polígono Bakiola, 41  
48498 Arrankudiaga  
Vizcaya - Spain

Tel: +34 (0) 902142434  
[www.quickiepower.com](http://www.quickiepower.com)  
For international enquiries e-mail:  
[International@sunrisemedical.de](mailto:International@sunrisemedical.de)

All information is subject to change without notification.  
Please consult Sunrise Medical with any queries you may have.  
If you are visually impaired, this document can be viewed in PDF  
format at [www.SunriseMedical.com](http://www.SunriseMedical.com) or alternatively is available on  
request in large text.